


How to install driver for AutoTiss One (Win 10)


1. Turn off your internet
2. Power off AutoTiss One
3. DO NOT plug-in AutoTiss One
4. Double click install_driver.exe to install driver
5. Plug-in AutoTiss One to your computer







6. Go to search icon at the lower left corner
7. Type in “device manager” to launch the driver update

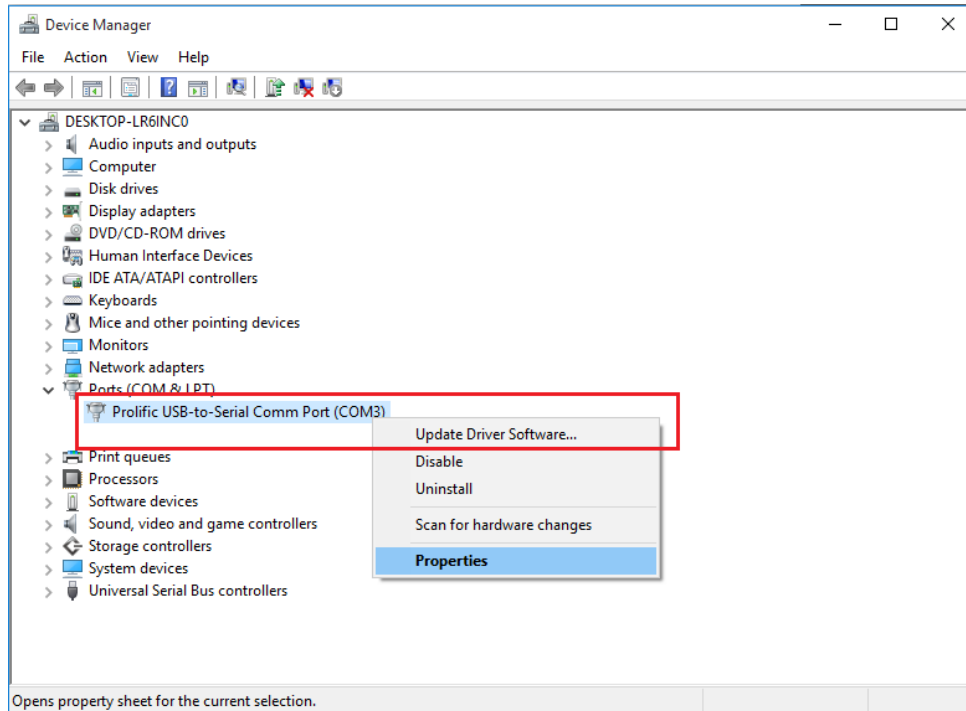

8. Right mouse button click on “Prolific USB-to-Serial Comm Port (COM)” and select “Update Driver Software...”

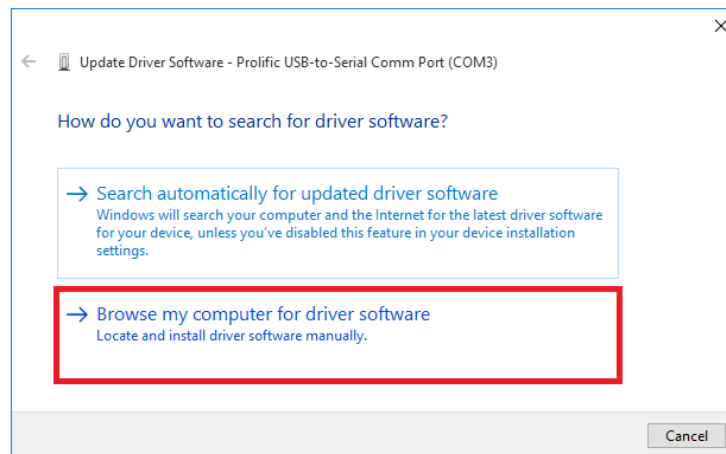

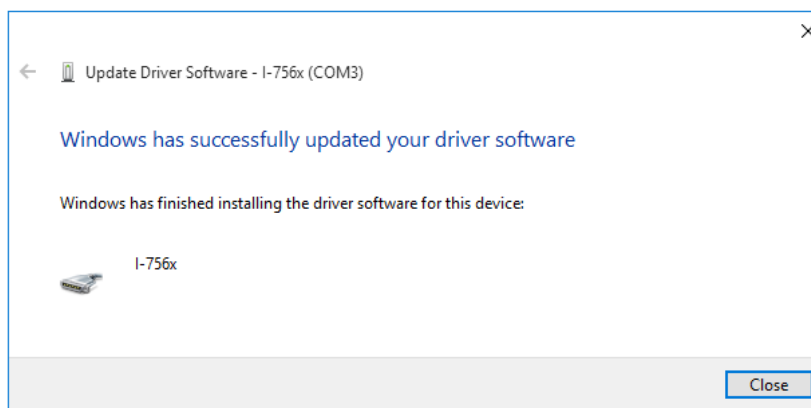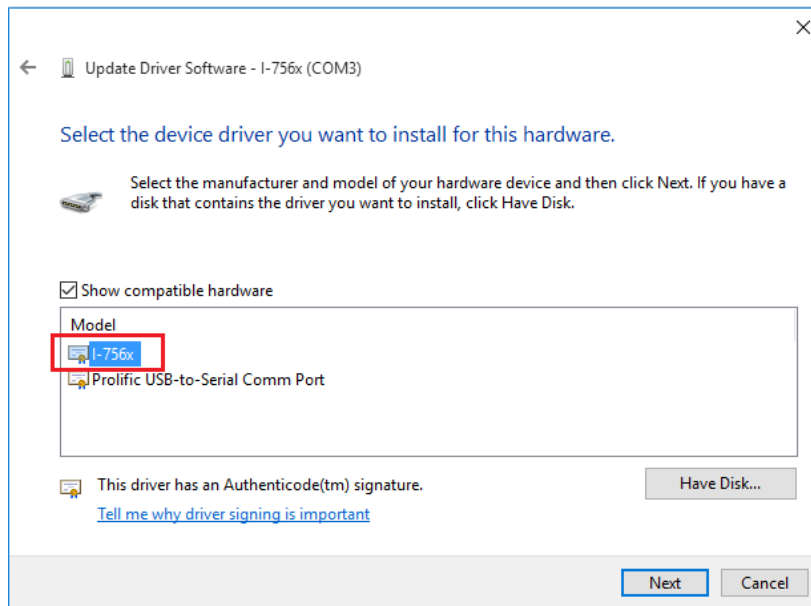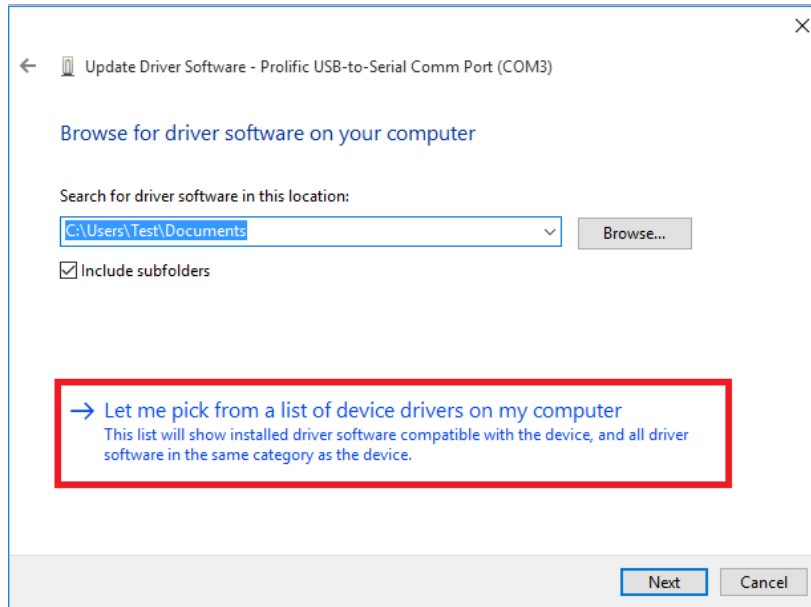

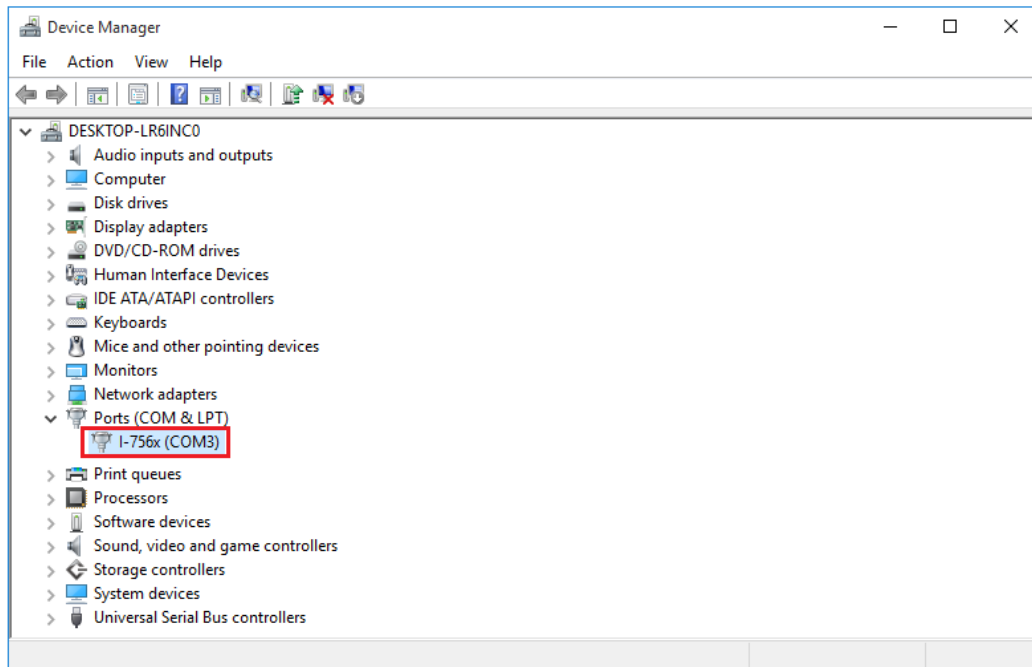

9. COM port number might be different on your computer. Open the setting.txt file and change accordingly.

```
com 3
dirX 0
dirY 0
homeX 0
homeY 0
homeDirX 0
homeDirY 0
ppmX 51200
ppmY 51200
mprX 2
mprY 2
homeOffsetX 2100000
homeOffsetY 700000
limitX 80
limitY 25
donorX 64
donorY 12
```