





How to install driver for AutoTiss One (Win 7)



1. Turn off your internet
2. Power off AutoTiss One
3. DO NOT plug-in AutoTiss One
4. Copy AutoTissOne folder in the USB disk to your computer desktop
5. Double click install_driver_win7.exe to install driver
6. Go to device manager” to launch the driver update

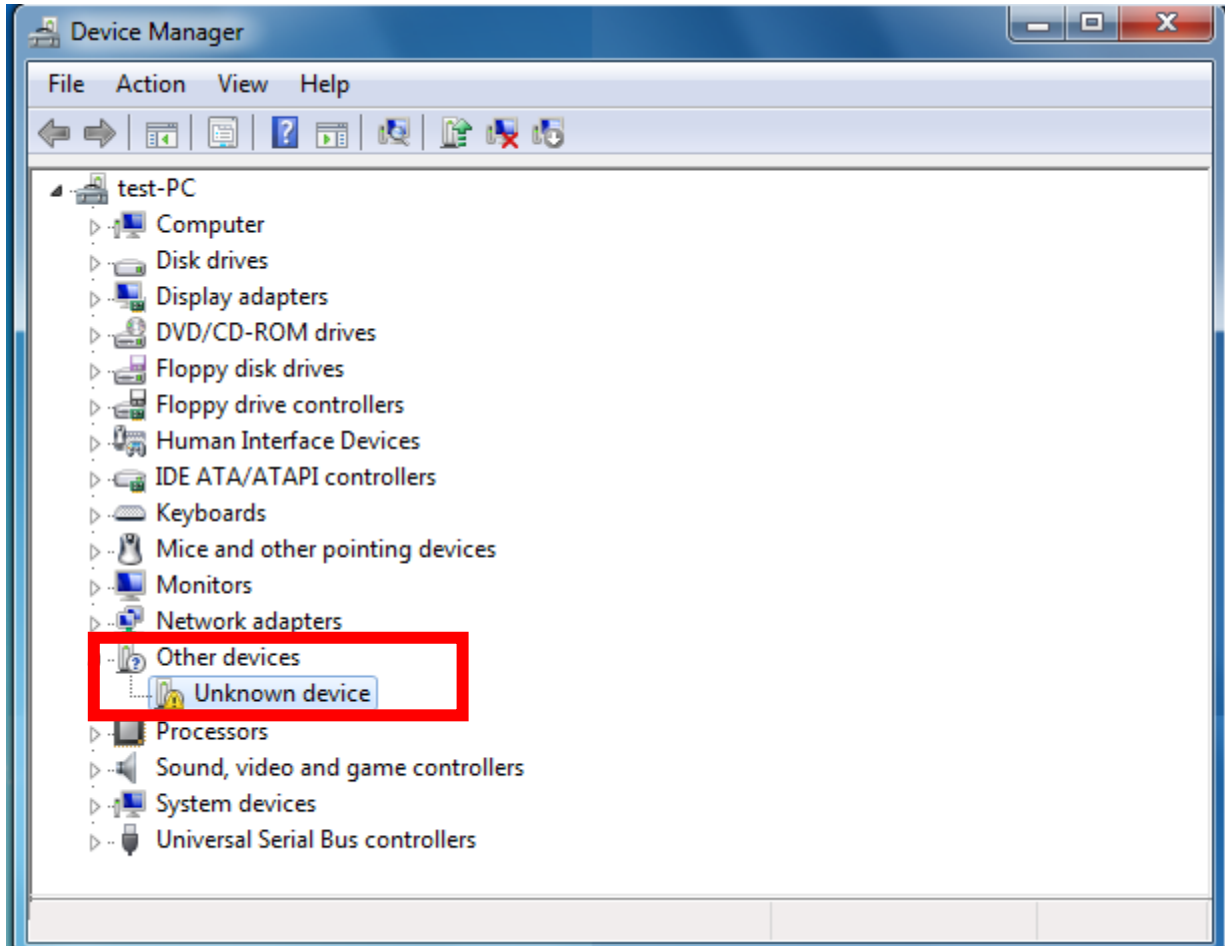


7. Right mouse button click on “Unknown device” and select “properties”

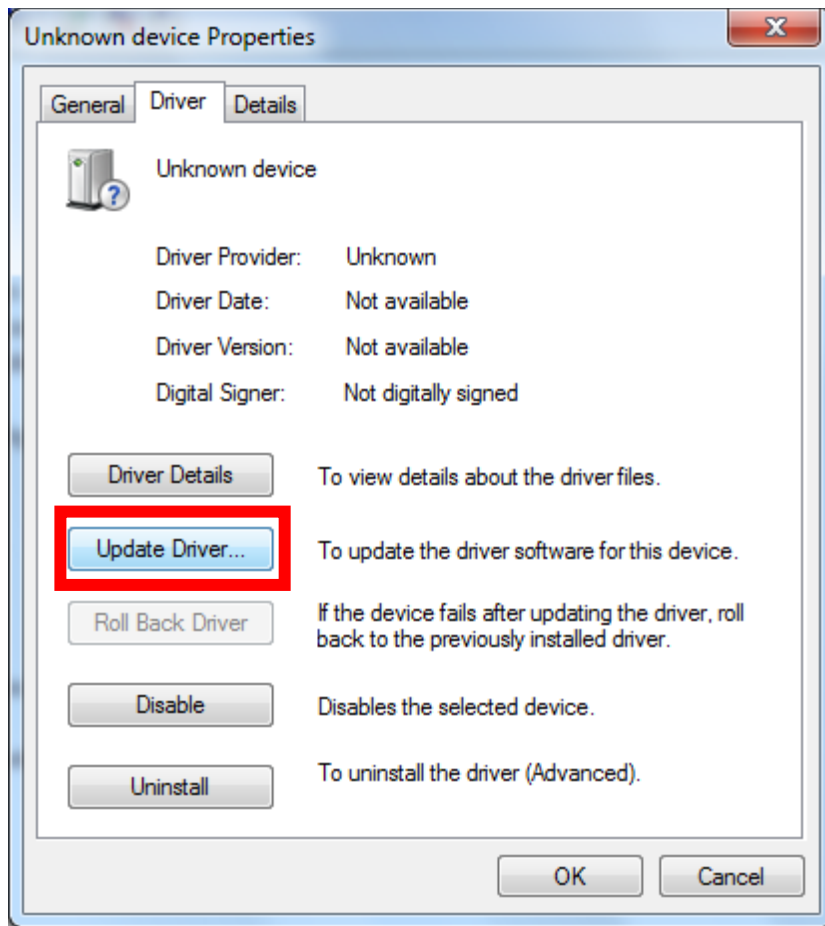

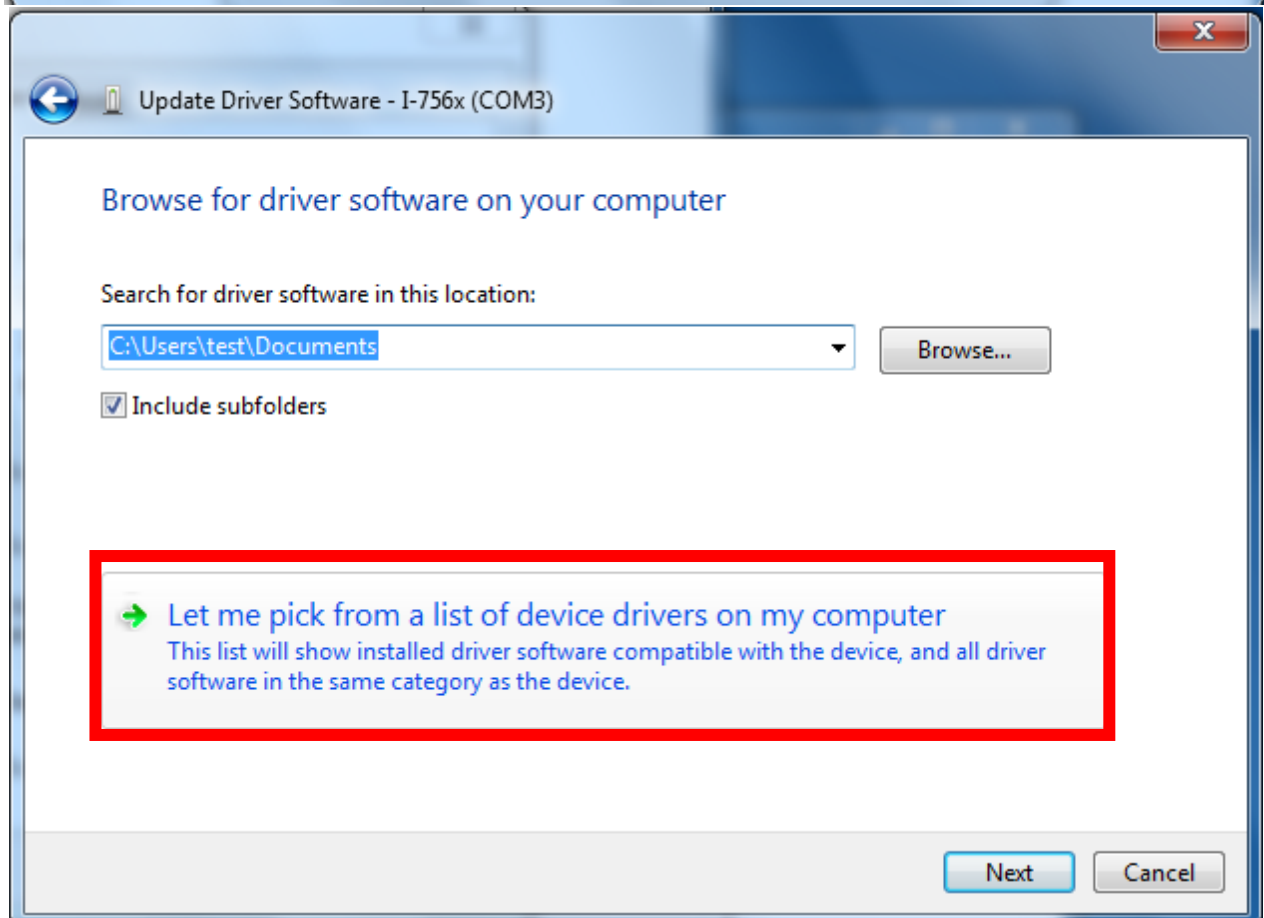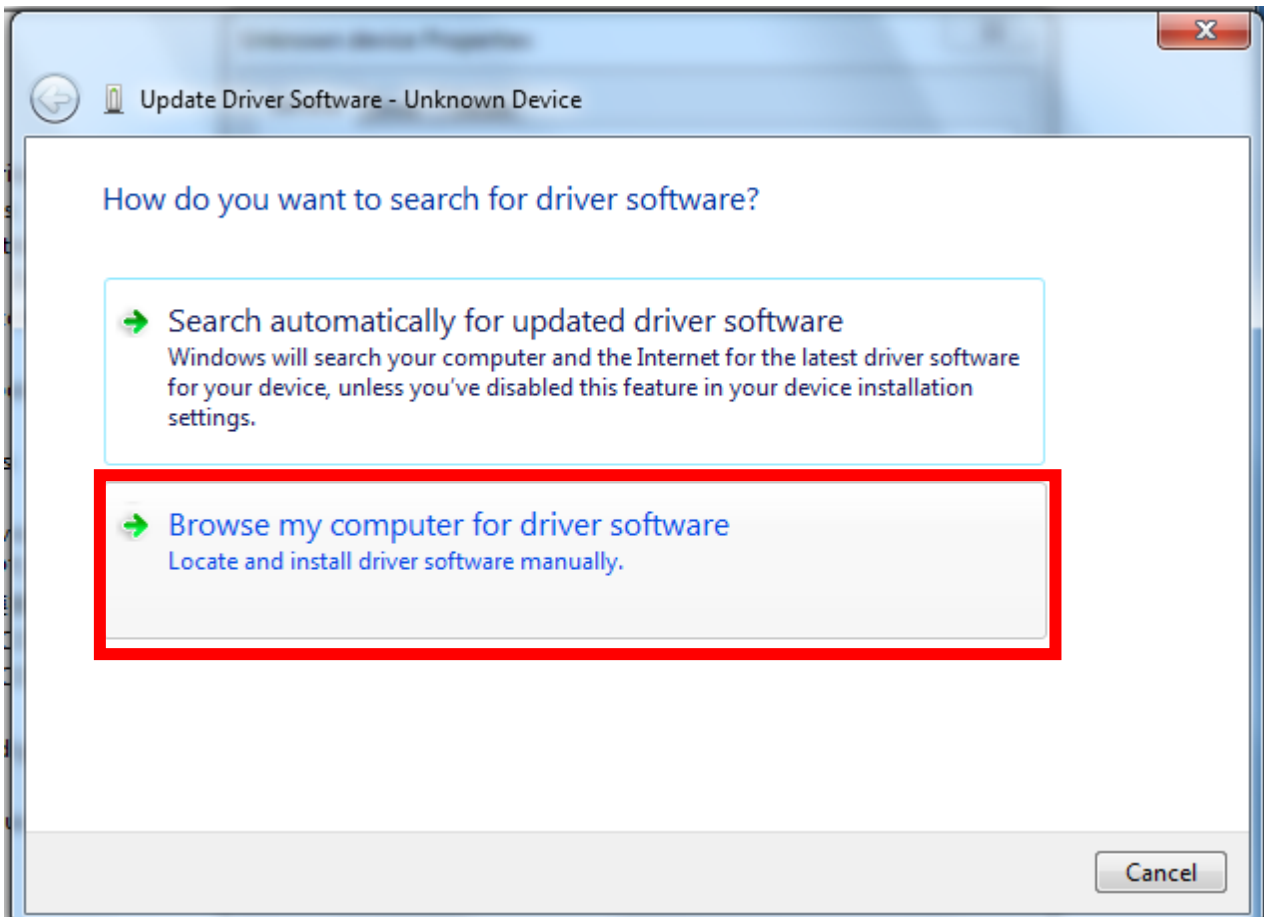

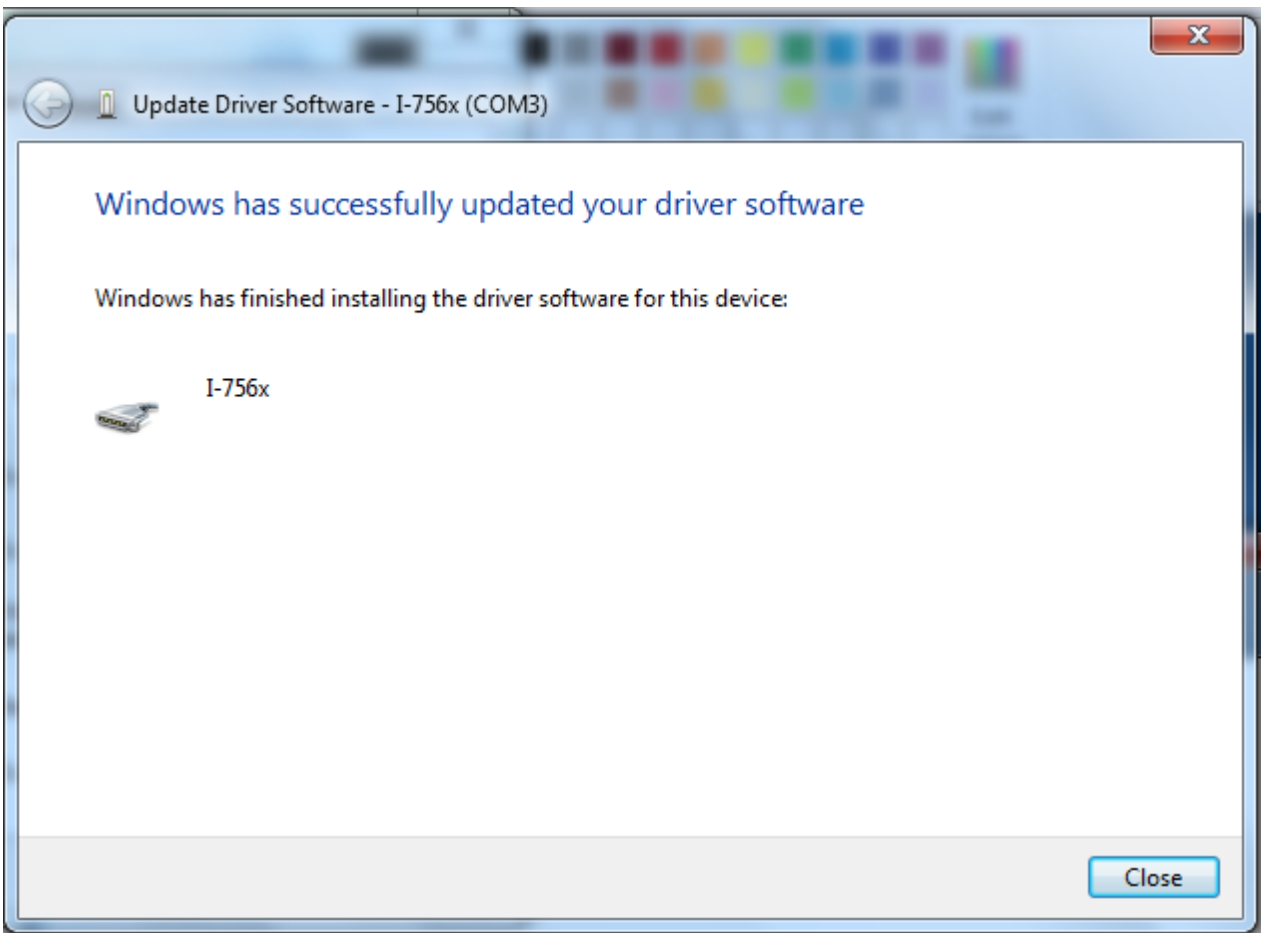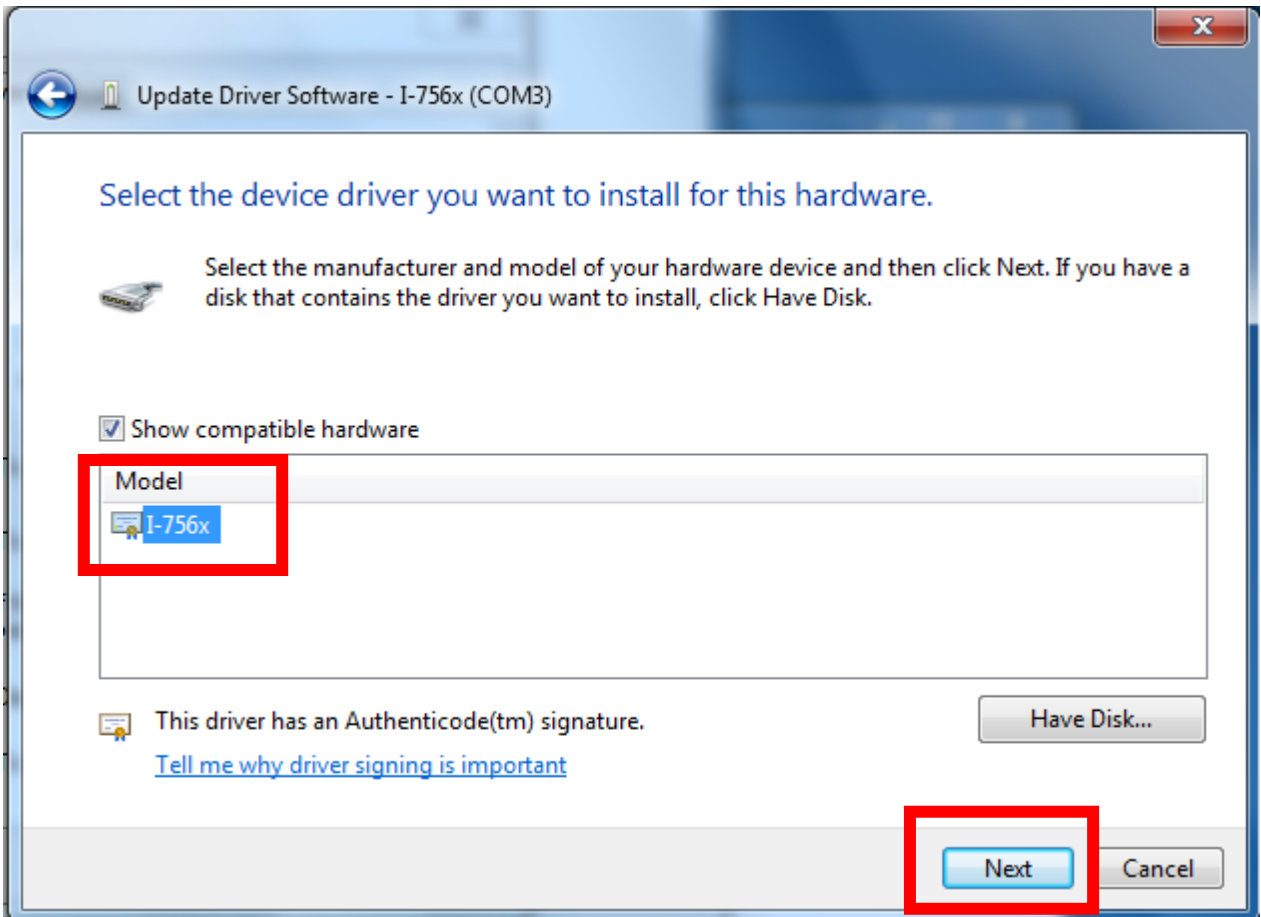

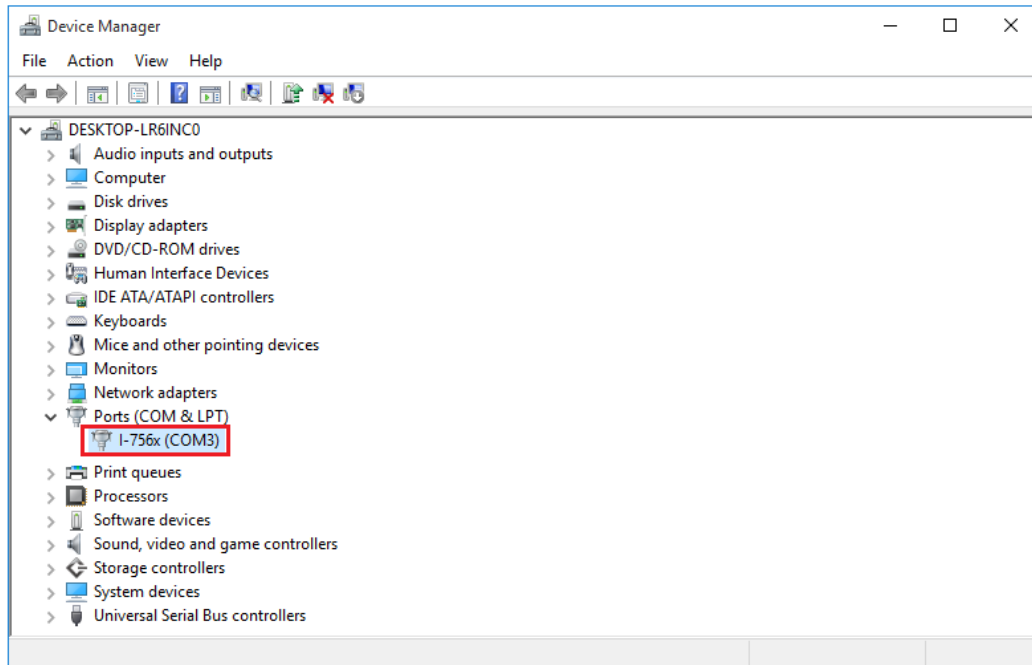

8. COM port number might be different on your computer. Open the setting.txt file and change accordingly.

```
com 3  
dirX 0  
dirY 0  
homeX 0  
homeY 0  
homeDirX 0  
homeDirY 0  
ppmX 51200  
ppmY 51200  
mprX 2  
mprY 2  
homeOffsetX 2100000  
homeOffsetY 700000  
limitX 80  
limitY 25  
donorX 64  
donorY 12
```

9. Save and close the setting file. You are now ready to start the program AutoTissOne.exe